

# icare

## IC333

The IC333 Bed is unmatched for function and durability, taking the hospital look out of your room. This four function bed, features dual lift motors totalling 12,000 newtons lifting power, unlimited under-bed clearance and eight heavy duty castors.



### KEY FUNCTIONS:



Head Lift



Knee Break



Height Adjustable  
22-66cm



Trendelenburg and  
Reverse Trendelenburg

### BED CODES & SIZES:



LONG SINGLE  
Code: IC333LS  
Mattress Surface Size: 2030 x 900  
Base Size: 2150 x 960



LONG DOUBLE  
Code: IC333D  
Mattress Surface Size: 2030 x 1350  
Base Size: 2150 x 1400  
One hand control



KING SINGLE  
Code: IC333KS  
Mattress Surface Size: 2030 x 1070  
Base Size: 2150 x 1120




QUEEN  
Code: IC333Q  
Mattress Surface Size: 2030 x 1500  
Base Size: 2150 x 1580  
One hand control

Note: Always remember that TWO Long Single beds together make a King size bed.

## KEY FEATURES:



200kg SWL on Long Single And King Single

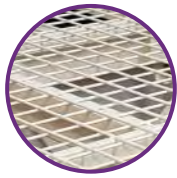

Strong Breathable Mattress Platform



Available With Headboards



Rechargeable Battery Backup



Lockable Heavy Duty Castors



300kg SWL on Double and Queen



End and Side Mattress Retainer Bars



10 Year Warranty



Back Lit Wipeable Hand Piece



No Gap Between Mattresses on King Beds (2x long singles)



Unlimited Under Bed Clearance



Ultra Low Design Reducing Risk of Falls

## TECHNICAL SPECIFICATIONS:

- » Safe working load - 200kg on long single and king single sizes and 300kg for double and queen sizes. The safe working load includes any model of Icare mattresses.
- » Hand Piece - Backlit for night time use. Easy press buttons. Wired stretch cord. Hook attachment for hanging to side of bed. Wipeable.
- » Fabric - Hard wearing commercial grade fabric. Onyx range (BLACK) 82% Polyester, 18% Olefin. Stone range (GREY) 70% Polyester, 30% Acrylic.
- » Mattress Surface - Strong metal mesh under mattress allowing increased airflow
- » Actuators - 1x 6000 newton - 24 watt head adjust actuator. 1x 6000 newton - 24 watt knee adjust actuator. 2x 6000 newton - 24 watt height adjust actuators (2x 8000 newton on double and queen sizes)
- » Wheels - 8x heavy duty swivel castors with one lockable castor on each corner.
- » Head lift angle 62° and knee break angle 42°
- » Max height - 66cm from floor to top of base
- » Min height - 22cm from floor to top of base

## ACCESSORIES AVAILABLE:

- » High Side Rail
- » Low Side Rail
- » Bed Stick
- » Over Bed Pole
- » Full Length Side Rails
- » Headboards



Available in one piece Double and Queen